

HPLC

Instrument: Waters-HPLC **Column:** C18, 4.6*100 mm


Detection Mode: 265nm **Column temperature:** 30°C

Flow Rate: 1.0ml/min **Sample dissolution:** Methanol

Mobile Phase: A : Acetonitrile B : 0.1%phosphoric acid in water;

Gradientelution: A: 45%-55%-55%, 0-15-25min

CAS NO.: 77029-83-5

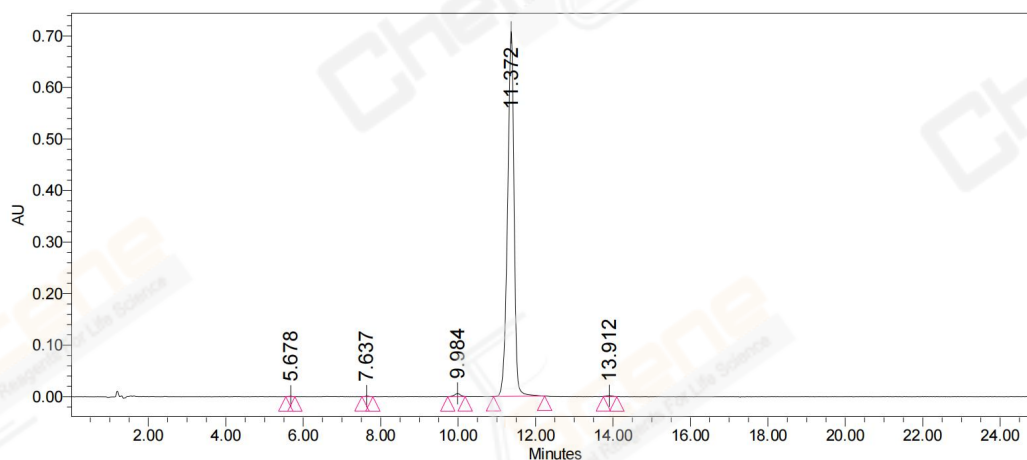

	Name	RT	Area	% Area	USP Plate Count	USP Resolution
1		5.678	6089	0.07	11266.8	
2		7.637	7641	0.09	8833.2	9.7
3		9.984	64286	0.77	21912.8	8.9
4		11.372	8301942	98.83	22544.3	4.5
5		13.912	20485	0.24	24605.6	7.0